

Certificate of Analysis for NR-3037

Genomic DNA from Francisella tularensis subsp. novicida, Strain CG116

Catalog No. NR-3037

Product Description: Genomic DNA was isolated from a preparation of *Francisella tularensis* (*F. tularensis*) subsp. *novicida*, strain CG116.

Lot¹: 7639520 Manufacturing Date: 29SEP2006

TEST	SPECIFICATIONS	RESULTS
Sequencing of 16S Ribosomal RNA Gene (~ 1380 bp)	Identical to GenBank CP000439 Consistent with F. tularensis subsp. novicida	Identical to GenBank CP000439 Consistent with <i>F. tularensis</i> subsp. <i>novicida</i> ²
Agarose Gel Electrophoresis	High molecular weight chromosomal DNA	High molecular weight chromosomal DNA (Figure 1)
Concentration by PicoGreen® Measurement	4 to 6 μg in 25 to 100 μL per vial	5.3 μg in 55 μL per vial (97 μg/mL)
Functional Activity by PCR Amplification 16S ribosomal RNA gene Molecular subtyping ^{3,4}	~ 1500 bp amplicon ~ 390 bp amplicon (non-type B)	~ 1500 bp amplicon ~ 390 bp amplicon (non-type B)
OD ₂₆₀ /OD ₂₈₀ Ratio	1.7 to 1.9	1.9
Bacterial Inactivation 10% of total yield plated on Trypticase Soy Agar with 5% sheep blood ^{5,6}	No viable bacteria detected	No viable bacteria detected

¹F. tularensis subsp. novicida, strain CG116 was extracted from NR-582 (Lot: 4067737) on Difco™ Cystine Heart Agar (BD 247100) with 5% defibrinated rabbit blood. After incubation for 24 hours at 37°C in an aerobic atmosphere with 5% CO₂, genomic DNA was extracted using proprietary technology.

Date: 24 JUL 2007 **Signature:** Signature on File

Title: Technical Manager, BEI Authentication or designee

ATCC[®], on behalf of BEI Resources, hereby represents and warrants that the material provided under this certificate has been subjected to the tests and procedures specified and that the results described, along with any other data provided in this certificate, are true and accurate to the best of ATCC[®]'s knowledge.

ATCC® is a trademark of the American Type Culture Collection.

You are authorized to use this product for research use only. It is not intended for human use.

Fax: 703-365-2898

E-mail: contact@beiresources.org

800-359-7370

Also consistent with other *F. tularensis* subspecies.

³Petersen, J. M., et al. "Laboratory Analysis of Tularemia in Wild-Trapped, Commercially Traded Prairie Dogs, Texas, 2002." <u>Emerg. Infect. Dis.</u> 10 (2004): 419–425. PubMed: 15109407.


⁴Kugeler, K. J., et al. "Real-time PCR for Francisella tularensis Types A and B." Emerg. Infect. Dis. 12 (2006): 1799–1801. PubMed: 17283646.

⁵Incubated for 8 days at 37°C and aerobic atmosphere with 5% CO₂.

⁶An extraction procedure was used that has been shown to consistently inactivate 100% of Gram-negative bacteria.



Certificate of Analysis for NR-3037

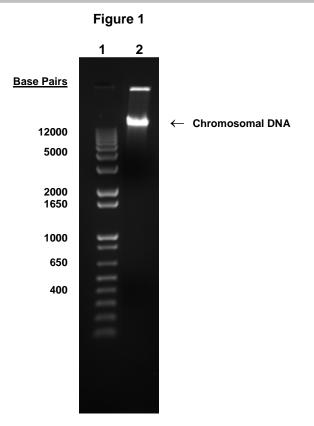

Lane 1: Invitrogen™ TrackIt™ 1 Kb Plus DNA Ladder

Lane 2: 200 ng of NR-3037

Fax: 703-365-2898

E-mail: contact@beiresources.org

800-359-7370