

SUPPORTING INFECTIOUS DISEASE RESEARCH

Certificate of Analysis for NR-10043

Genomic RNA from Influenza A Virus, A/New Caledonia/20/99 (HA, NA) x A/Puerto Rico/8/34 (H1N1), Reassortant X-139 (Kilbourne F98)

Catalog No. NR-10043

Product Description: Genomic RNA was isolated from a preparation of pooled allantoic fluid from specific-pathogen free (SPF) embryonated chicken eggs¹ infected with reassortant influenza A virus, A/New Caledonia/20/99 (HA, NA) x A/Puerto Rico/8/34 (H1N1) (Kilbourne F98; X-139).

Lot²: 58330565 Manufacturing Date: 12JUN2009

TEST	SPECIFICATIONS	RESULTS
Sequencing of Influenza A Specific Region Matrix gene (776 nucleotides)	Identical to BEI Resources NR-9692 Consistent with influenza A virus	Identical to BEI Resources NR-9692 Consistent with influenza A virus
Total RNA Content by RiboGreen® Measurement (Viral, Cellular and Carrier)	Report results	9.1 ng per 100 µL
Functional Activity by RT-PCR Amplification ^{3,4}	~ 1030 bp amplicon	~ 1030 bp amplicon (See Figure 1)
Virus Inactivation 10% of total yield incubated in SPF embryonated chicken eggs ¹ and evaluated by hemagglutination assay ⁵	No virus detected	No virus detected
Sodium Azide Content	Report results	≤ 0.0004%

¹⁰ to 11-day-old SPF Embryonated Chicken Eggs acquired from B&E Eggs, York Springs, Pennsylvania

Date: 18 JAN 2012 **Signature:**

Title: Technical Manager, BEI Authentication or designee

Dorothy C. Young

ATCC[®], on behalf of BEI Resources, hereby represents and warrants that the material provided under this certificate has been subjected to the tests and procedures specified and that the results described, along with any other data provided in this certificate, are true and accurate to the best of ATCC[®]'s knowledge.

ATCC® is a trademark of the American Type Culture Collection.

You are authorized to use this product for research use only. It is not intended for human use.

BEI Resources www.beiresources.org contact@beiresources.org
Tel: 800-359-7370

Fax: 703-365-2898
NR-10043_58330565_18JAN2012

²Nucleic acid was extracted from a preparation of reassortant influenza A virus, A/New Caledonia/20/99 (HA, NA) x A/Puerto Rico/8/34 (H1N1), Reassortant X-139 (Kilbourne F98) (BEI Resources NR-9692; Lot 58155270) using a QIAamp[®] Viral RNA Mini kit (Qiagen 52906).

³Amplified using 1:10 dilution of NR-10043 and One-Step RT-PCR Kit (Qiagen 210212)


⁴The primers are described in Hoffmann, E., et al. "Universal Primer Set for the Full-Length Amplification of All Influenza A Viruses." <u>Arch. Virol.</u> 146 (2001): 2275-2289. PubMed: 11811679.

⁵This extraction procedure has been shown to consistently inactivate 100% of influenza viruses.



SUPPORTING INFECTIOUS DISEASE RESEARCH

Certificate of Analysis for NR-10043

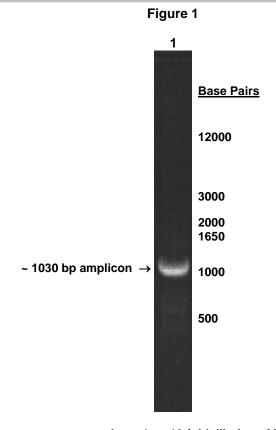

Lane 1: 10-fold dilution of NR-10043

Markers: Invitrogen™ TrackIt™ 1 Kb Plus DNA Ladder