

Certificate of Analysis for NR-13601

Genomic DNA from Clostridium difficile, Isolate 10

Catalog No. NR-13601

Product Description: Genomic DNA was obtained from a preparation of *Clostridium difficile* (*C. difficile*), isolate 10 which was isolated from a human patient from the Mid-Atlantic region of the United States in 2008/2009.

Lot¹: 59147496 Manufacturing Date: 24AUG2010

TEST	SPECIFICATIONS	RESULTS
Sequencing of 16S Ribosomal RNA Gene (~ 1330 base pairs)	Consistent with C. difficile	Consistent with C. difficile
Agarose Gel Electrophoresis	High molecular weight chromosomal DNA	High molecular weight chromosomal DNA (Figure 1)
Concentration by PicoGreen® Measurement	0.7 to 1.5 μg in 25 to 100 μL per vial	0.8 μg in 40 μL per vial (21 μg/mL)
Functional Activity by PCR Amplification 16S ribosomal RNA gene	~ 1500 bp amplicon	~ 1500 bp amplicon
OD ₂₆₀ /OD ₂₈₀ Ratio	1.7 to 1.9	1.7
Bacterial Inactivation 10% of total yield plated on CDC anaerobic blood agar ^{2,3}	No viable bacteria detected	No viable bacteria detected

¹The bacterial preparation used for extraction of genomic DNA was produced by Modified Reinforced Clostridial Broth (<u>ATCC medium 2107</u>) culture of the deposited material. After incubation for 48 hours at 37°C in an anaerobic atmosphere (80% N₂:10% CO₂:10% H₂), genomic DNA was extracted using proprietary technology.

Date: 11 FEB 2011 Signature:

Title: Technical Manager, BEI Authentication or designee

ATCC®, on behalf of BEI Resources, hereby represents and warrants that the material provided under this certificate has been subjected to the tests and procedures specified and that the results described, along with any other data provided in this certificate, are true and accurate to the best of ATCC®'s knowledge.

ATCC® is a trademark of the American Type Culture Collection.

You are authorized to use this product for research use only. It is not intended for human use.

Biodefense and Emerging Infections Research Resources Repository

www.beiresources.org

Tel: 800-359-7370


NR-13601_59147496_11FEB2011

²7 days at 37°C in an anaerobic atmosphere (80% N₂:10% CO₂:10% H₂)

³An extraction procedure was used that has been shown to consistently inactivate 100% of Gram-negative and Gram-positive bacteria.



Certificate of Analysis for NR-13601

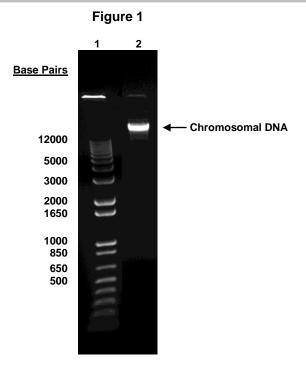

Lane 1: Invitrogen™ TrackIt 1 Kb Plus DNA Ladder™

Lane 2: 200 ng of NR-13601

NR-13601_59147496_11FEB2011