

Certificate of Analysis for NR-46024

Staphylococcus aureus, Strain HT 20020058

Catalog No. NR-46024

Product Description: *Staphylococcus aureus* (*S. aureus*), strain HT 20020058 was isolated in 2002 from the arthritic knee joint of a 7-year-old female in France. *S. aureus*, strain HT 20020058 is a clinically associated methicillin-sensitive *S. aureus* (MSSA) strain.

Lot¹: 63622050 Manufacturing Date: 10JUL2015

TEST	SPECIFICATIONS	RESULTS
Phenotypic Analysis		
Cellular morphology	Gram-positive cocci	Gram-positive cocci
Colony morphology ²	Report results	Circular, flat, entire, smooth and cream
, , ,	·	(Figure 1)
Motility (wet mount)	Report results	Non-motile
Hemolysis ²	Report results	β-hemolytic
Biochemical Analysis		
Catalase	Positive	Positive
Coagulase ³	Report results	Positive
VITEK [®] 2 Compact (GP card)	Consistent with S. aureus	Consistent with S. aureus
Antibiotic Susceptibility Profile		
VITEK [®] (AST-GP71 card) ⁴		
Beta-lactamase ⁵	Report results	Positive
Cefoxitin screen	Report results	Negative
Benzylpenicillin	Report results	Resistant (≥ 0.5 μg/mL)
Oxacillin	Sensitive	Sensitive (≤ 0.5 μg/mL)
Gentamicin	Sensitive	Sensitive (≤ 0.5 μg/mL)
Ciprofloxacin	Sensitive	Sensitive (≤ 0.5 µg/mL)
Levofloxin	Report results	Sensitive (≤ 0.12 μg/mL)
Moxifloxacin	Report results	Sensitive (≤ 0.25 μg/mL)
Clindamycin (inducible resistance)	Report results	Negative
Erythromycin	Sensitive	Sensitive (≤ 0.25 μg/mL)
Clindamycin	Sensitive	Sensitive (≤ 0.25 μg/mL)
Quinupristin/dalfopristin	Sensitive	Sensitive (≤ 0.25 μg/mL)
Linezolid	Sensitive	Sensitive (= 2 µg/mL)
Daptomycin	Report results	Sensitive (≤ 0.5 μg/mL)
Vancomycin Minocycline	Sensitive Report results	Sensitive (≤ 0.5 µg/mL) Sensitive (≤ 0.5 µg/mL)
Tetracycline	Report results	Sensitive (≤ 0.5 μg/mL)
Tigecycline	Report results	Sensitive ($\leq 1 \mu g/mL$) Sensitive ($\leq 0.12 \mu g/mL$) ⁶
Nitrofurantoin	Report results	Sensitive (≤ 32 µg/mL)
Rifampicin	Report results	Sensitive (≤ 0.5 µg/mL)
Trimethoprim/sulfamethoxazole	Sensitive	Sensitive (≤ 0.5 µg/mL)
Etest [®] antibiotic test strips ⁷		(- · · · · · · · · · · · · · · · ·
Chloramphenicol ⁸	Report results	Sensitive (= 3 µg/mL)
Teicoplanin ⁸	Sensitive	Sensitive (= 1 µg/mL)
Genotypic Analysis		
Sequencing of 16S ribosomal RNA gene	Consistent with S. aureus	Consistent with S. aureus
(~ 1500 base pairs)		
Purity (post-freeze) ⁹	Consistent with expected colony	Consistent with expected colony
	morphology	morphology
Viability (post-freeze) ²	Growth	Growth
	I and the second	I .

BEI Resources www.beiresources.org E-mail: contact@beiresources.org

Tel: 800-359-7370 Fax: 703-365-2898




SUPPORTING INFECTIOUS DISEASE RESEARCH

Certificate of Analysis for NR-46024

¹S. aureus, strain HT 20020058 was deposited to BEI Resources as part of the NARSA collection. NR-46024 was produced by inoculation of the deposited material into Tryptic Soy broth and grown 1 day at 37°C in an aerobic atmosphere. Broth inoculum was added to Tryptic Soy agar with 5% defibrinated sheep blood kolles which were grown 1 day at 37°C in an aerobic atmosphere to produce this lot.

²1 day at 37°C in an aerobic atmosphere on Tryptic Soy agar with 5% defibrinated sheep blood

 $^{^9}$ Purity of this lot was assessed for $\bar{7}$ days on Tryptic Soy agar with 5% defibrinated sheep blood at 37°C in an aerobic atmosphere.

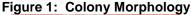



Date: 16 NOV 2015

Signature:

BEI Resources Authentication

ATCC®, on behalf of BEI Resources, hereby represents and warrants that the material provided under this certificate has been subjected to the tests and procedures specified and that the results described, along with any other data provided in this certificate, are true and accurate to the best of ATCC®'s knowledge.

ATCC® is a trademark of the American Type Culture Collection.

You are authorized to use this product for research use only. It is not intended for human use.

BEI Resources www.beiresources.org E-mail: contact@beiresources.org
Tel: 800-359-7370

Fax: 703-365-2898

³4 hours at 37°C in rabbit serum with 0.15% EDTA (Coagulase Plasma BBL™ 240827)

⁴Minimum Inhibitory Concentration (MIC); MIC Interpretation Guideline: CLSI M100-S22 (2012)

⁵The production of beta-lactamase was detected using a Cefinase™ Paper Disc (BBL™ 231650).

⁶MIC Interpretation Guideline: EUCAST Version 4.0 (2014)

⁷¹ day at 37°C in an aerobic atmosphere on Mueller Hinton agar

⁸For both chloramphenicol (bioMérieux Etest[®] 412308) and teicoplanin (bioMérieux Etest[®] 412459), a MIC ≤ 8 μg/mL is sensitive, a MIC = 16 μg/mL is intermediate and a MIC ≥ 32 μg/mL is resistant.