

Certificate of Analysis for NR-4612

Genomic DNA from *Salmonella enterica* subsp. *enterica*, 2004 Pennsylvania Tomato Outbreak, Serovar Thompson, Isolate 13

Catalog No. NR-4612

Product Description: Genomic DNA was obtained from a preparation of *Salmonella enterica* (*S. enterica*) subsp. *enterica* serovar Thompson that was isolated from the stool of a patient with diarrhea during the 2004 *Salmonella* outbreak in Pennsylvania.

Lot¹: 58332161 Manufacturing Date: 08OCT2008

TEST	SPECIFICATIONS	RESULTS
Sequencing of 16S Ribosomal RNA Gene (~ 1420 bp)	Identical to BEI Resources NR-4331 Consistent with <i>S. enterica</i>	Identical to BEI Resources NR-4331 Consistent with <i>S. enterica</i>
Agarose Gel Electrophoresis	High molecular weight chromosomal DNA	High molecular weight chromosomal DNA (Figure 1)
Concentration by PicoGreen® Measurement	4 to 6 μg in 25 to 100 μL per vial	5.4 μg in 38 μL per vial (141 μg/mL)
Functional Activity by PCR Amplification 16S ribosomal RNA gene	~ 1500 bp amplicon	~ 1500 bp amplicon
OD ₂₆₀ /OD ₂₈₀ Ratio	1.7 to 1.9	1.8
Bacterial Inactivation 10% of total yield plated on Tryptic Soy Agar ^{2,3}	No viable bacteria detected	No viable bacteria detected

¹S. enterica subsp. enterica, 2004 Pennsylvania Tomato Outbreak, Serovar Thompson, Isolate 13 was deposited by Carol H. Sandt, Molecular Microbiology Section, Bureau of Laboratories, PA Department of Health, Lionville, PA. The bacterial preparation used for extraction of genomic DNA was produced by broth (Tryptic Soy Broth, BD 211825) culture of the deposited material. After incubation for 24 hours at 37°C and aerobic atmosphere, genomic DNA was extracted using proprietary technology.

Date: 12 DEC 2008 **Signature:** Signature on File

Title: Technical Manager, BEI Authentication or designee

ATCC®, on behalf of BEI Resources, hereby represents and warrants that the material provided under this certificate has been subjected to the tests and procedures specified and that the results described, along with any other data provided in this certificate, are true and accurate to the best of ATCC®'s knowledge.

© 2008 American Type Culture Collection (ATCC). All rights reserved.

Page 1 of 2

 $\mbox{ATCC}^{\circledast}$ is a trademark of the American Type Culture Collection.

You are authorized to use this product for research use only. It is not intended for human use.

Biodefense and Emerging Infections Research Resources Repository P.O. Box 4137

Manassas, VA 20108-4137 USA www.beiresources.org

800-359-7370

Fax: 703-365-2898


NR-4612_58332161_12DEC2008

²7 days at 37°C in an aerobic atmosphere.

³An extraction procedure was used that has been shown to consistently inactivate 100% of Gram-negative bacteria.



Certificate of Analysis for NR-4612

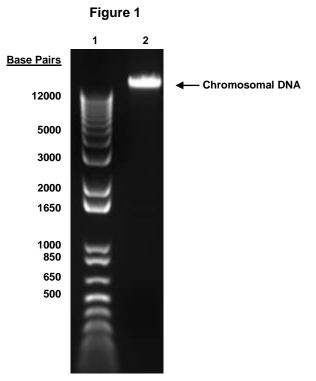

Lane 1: Invitrogen™ TrackIt 1 Kb Plus DNA Ladder™

Lane 2: 200 ng of NR-4612

Fax: 703-365-2898

E-mail: contact@beiresources.org

800-359-7370