

***Staphylococcus aureus* subsp. *aureus*, Strain JE2, Transposon Mutant NE27 (SAUSA300_1895)**

Catalog No. NR-46570

Product Description: *Staphylococcus aureus* (*S. aureus*) subsp. *aureus*, transposon mutant NE27 was derived from *S. aureus* subsp. *aureus*, strain JE2. Mutagenesis occurred through the use of the *mariner*-based transposon *bursa aurealis* resulting in an erythromycin-resistant deletion strain of JE2. *S. aureus* subsp. *aureus*, transposon mutant NE27 was created by disruption of SAUSA300_1895, which encodes for a nitric oxide synthase oxygenase that catalyzes nitric oxide production from L-arginine and may contribute to methicillin-resistant *S. aureus* (MRSA) innate immune and antibiotic resistance phenotypes. Strain JE2 is a plasmid-cured derivative of strain LAC that was isolated in 2002 from a skin and soft tissue infection of an inmate in the Los Angeles County Jail in California, USA.

Lot¹: 63303084

Manufacturing Date: 11FEB2015

TEST	SPECIFICATIONS	RESULTS
Phenotypic Analysis Cellular morphology Colony morphology ² Motility (wet mount)	Gram-positive cocci Report results Report results	Gram-positive cocci Circular, convex, entire, smooth and yellow (Figure 1) Non-motile
Confirmation of Transposon Insertion³	Resistant to erythromycin	Resistant to erythromycin
Purity (post-freeze)⁴	Growth consistent with <i>S. aureus</i>	Growth consistent with <i>S. aureus</i>
Viability (post-freeze)²	Growth	Growth


¹NR-46570 was produced by inoculation of the deposited material into Tryptic Soy broth with 5 µg/mL erythromycin and incubated for 22 hours at 37°C in an aerobic atmosphere. Broth inoculum was added to Tryptic Soy agar with 5 µg/mL erythromycin kolles which were grown 20 hours at 37°C in an aerobic atmosphere to produce this lot.

²20 hours at 37°C in an aerobic atmosphere on Tryptic Soy agar with 5 µg/mL erythromycin

³Prior to initiating work, it is recommended that the presence and location of the transposon is confirmed. Gene specific primers should be paired with either the "Upstream" primer (5'-CTCGATTCTATTAACAAGGG-3') for transposons in the "plus" orientation or the "Buster" primer (5'-GCTTTTTCTAAATGTTTTTAAGTAAATCAAGTAC-3') for transposons in the "minus" orientation. For additional information, refer to Fey, P. D., et al. "A Genetic Resource for Rapid and Comprehensive Phenotype Screening of Nonessential *Staphylococcus aureus* Genes." *MBio* 4 (2013): e00537-12. PubMed: 23404398.

⁴Purity of this lot was assessed for 7 days at 37°C in an aerobic atmosphere on Tryptic Soy agar with 5% defibrinated sheep blood.

Figure 1

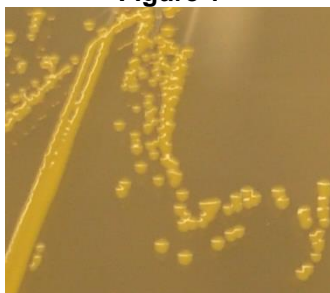

Date: 01 APR 2015

Signature: 

Title: Technical Manager, BEI Authentication or designee

ATCC[®], on behalf of BEI Resources, hereby represents and warrants that the material provided under this certificate has been subjected to the tests and procedures specified and that the results described, along with any other data provided in this certificate, are true and accurate to the best of ATCC[®]'s knowledge.

ATCC[®] is a trademark of the American Type Culture Collection.

You are authorized to use this product for research use only. It is not intended for human use.

