

Certificate of Analysis for NR-17787

Acinetobacter baumannii, Strain IS-123

Catalog No. NR-17787


Product Description: Acinetobacter baumannii (A. baumannii), strain IS-123 is a human isolate collected in February 2009 from the wound of a patient at Ibn Sina, Iraq.

Lot¹: 70018988 Manufacturing Date: 19SEP2018

TEST	SPECIFICATIONS	RESULTS
Phenotypic Analysis		
Cellular morphology	Gram-negative rods	Gram-negative rods
Colony morphologies ^{2,3}	Report results	Colony type 1: Circular, convex, entire, smooth and pale gray Colony type 2: Circular, low convex, entire, smooth and gray
Growth at 44°C ± 2°C ^{4,5}	Growth	Growth
Motility (wet mount)	Report results	Non-motile
Biochemical tests: Catalase	Positive	Positive
Oxidase	Negative	Negative
VITEK® Mass Spectrometry (MALDI-TOF)	Consistent with A. baumannii	Consistent with A. baumannii complex ⁶
Genotypic Analysis Sequencing of 16S ribosomal RNA gene (~ 1490 base pairs)	≥ 99% sequence identity to A. baumannii, strain IS-123 (GenBank: ALII01000002.1)	99.9% sequence identity to A. baumannii, strain IS-123 (GenBank: ALII01000002.1)
Purity (post-freeze) ⁷	Growth consistent with expected colony morphology	Growth consistent with expected colony morphology
Viability (post-freeze) ²	Growth	Growth

NR-17787 lot 70018988 was produced by inoculation of BEI Resources NRS-17787 lot 63305791 into Tryptic Soy broth and incubated for 1 day at 37°C in an aerobic atmosphere. Broth inoculum was added to Tryptic Soy agar kolles and incubated for 1 day at 37°C in an aerobic atmosphere to produce this lot.

Figure 1: Colony Morphology

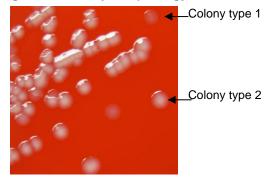

BEI Resources

www.beiresources.org

E-mail: contact@beiresources.org

Tel: 800-359-7370 Fax: 703-365-2898

²1 day at 37°C in an aerobic atmosphere on Tryptic Soy agar with 5% defibrinated sheep blood

³Two colony types were observed. Plating of the individual colony types showed that they did not revert to the mixed colony type. VITEK® MS (MALDITOF) analysis identified cells from both colony types as *A. baumannii*. The 16S ribosomal RNA gene of each colony type was sequenced and found to have 100% sequence identity to the other colony type and ≥ 99.9% sequence identity to *A. baumannii*, strain IS-123 (GenBank: ALII01000002.1)

⁴1 day at 44°C ± 2°C in an aerobic atmosphere on Tryptic Soy agar with 5% defibrinated sheep blood

⁵Growth at 44°C ± 2°C differentiates A. baumannii from A. calcoaceticus and A. pittii, which do not grow at 44°C ± 2°C.

⁶A. baumannii complex species include A. baumannii, A. calcoaceticus, A. pittii and A. nosocomialis.

⁷Purity of this lot was assessed for 7 days at 37°C in an aerobic atmosphere on Tryptic Soy agar with 5% defibrinated sheep blood.



Certificate of Analysis for NR-17787

/Heather Couch/ Heather Couch

27 DEC 2018

Program Manager or designee, ATCC Federal Solutions

ATCC®, on behalf of BEI Resources, hereby represents and warrants that the material provided under this certificate has been subjected to the tests and procedures specified and that the results described, along with any other data provided in this certificate, are true and accurate to the best of ATCC®'s knowledge.

ATCC® is a trademark of the American Type Culture Collection.

You are authorized to use this product for research use only. It is not intended for human use.

BEI Resources www.beiresources.org E-mail: contact@beiresources.org

Tel: 800-359-7370 Fax: 703-365-2898