


**Summary of Results**


Figure No.	Testing Conditions			Results
	Antibiotic Selection	Staining	Type of Stain	% of Cells Positive for CD4
1	Not selected	Not Stained	N/A	0.10%
2	Not selected	Isotype Stained	PE	0.02%
3	Not selected	CD4 Stained	PE	96.31%
4	Antibiotic selected	Not Stained	N/A	0.02%
5	Antibiotic selected	Isotype Stained	PE	0.05%
6	Not selected	CD4 Stained	PE	94.24%
7	Not selected	Isotype Stained	APC	0.03%
8	Not selected	CD4 Stained	APC	98.02%
9	Antibiotic selected	Isotype Stained	APC	0.06%
10	Antibiotic selected	CD4 Stained	APC	97.85%

Figure 1: HeLa CD4+ Cells (Clone 1022) were not selected with antibiotics and not stained

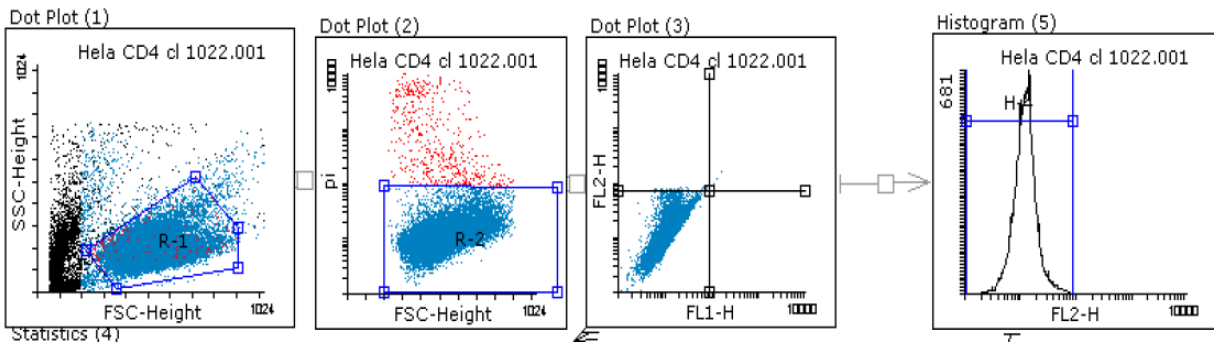


Results: 0.10% of cells showed positive fluorescence for CD4.

Figure 2: HeLa CD4+ Cells (Clone 1022) were not selected with antibiotics and isotype was PE stained (Isotype Control)

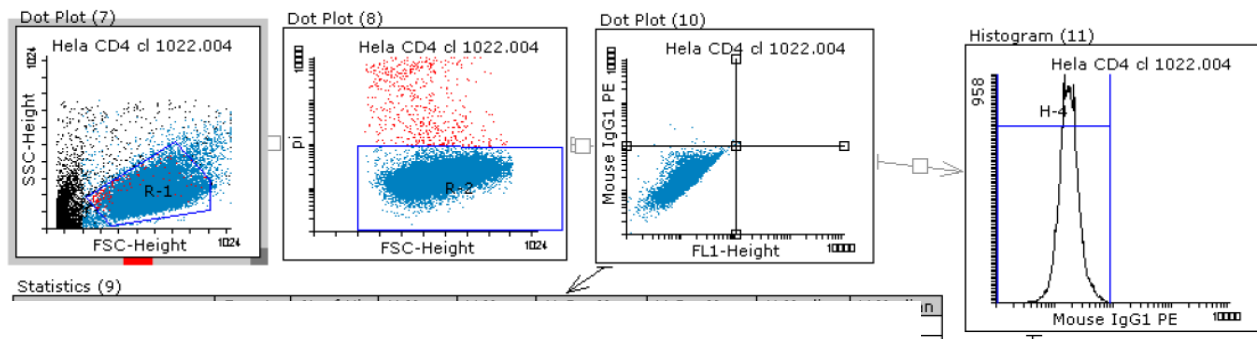


Results: 0.02% of cells showed positive fluorescence for CD4.

Figure 3: HeLa CD4+ Cells (Clone 1022) were not selected with antibiotics and CD4 was PE stained

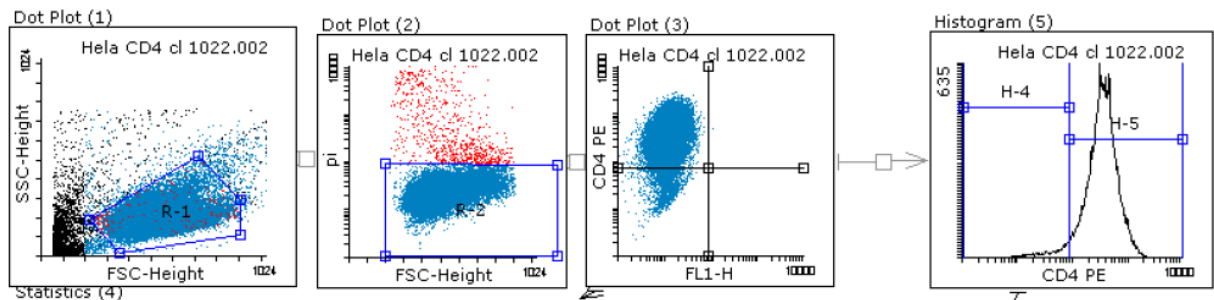


Results: 96.31% of cells showed positive fluorescence for CD4.

Figure 4: HeLa CD4+ Cells (Clone 1022) were selected with antibiotics and was not stained (Negative Control)

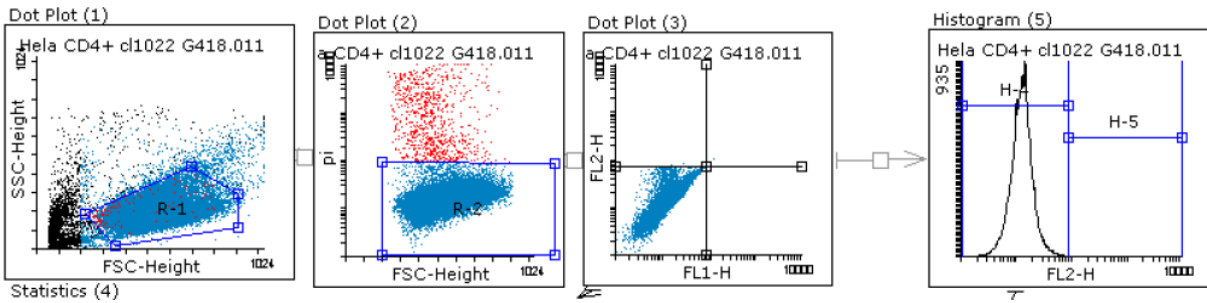


Results: 0.02% of cells showed positive fluorescence for CD4.

Figure 5: HeLa CD4+ Cells (Clone 1022) were selected with antibiotics and isotype was PE stained (Isotype Control)

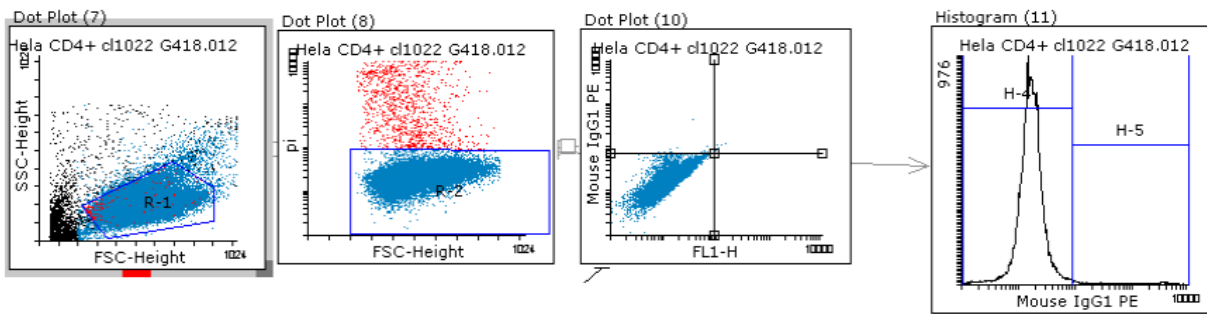


Results: 0.05% of cells showed positive fluorescence for CD4.

Figure 6: HeLa CD4+ Cells (Clone 1022) were selected with antibiotics and CD4 was PE stained.

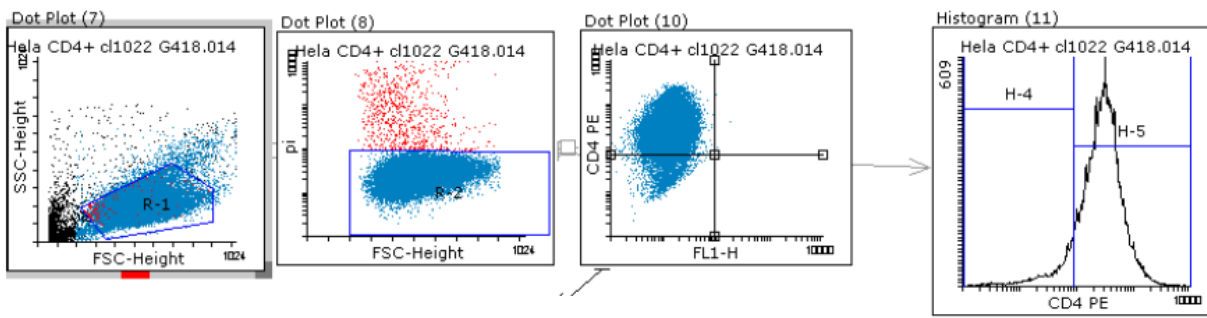


Results: 94.24% of cells showed positive fluorescence for CD4.

Figure 7: HeLa CD4+ Cells (Clone 1022) were not selected with antibiotics and isotype was APC stained (Isotype Control)

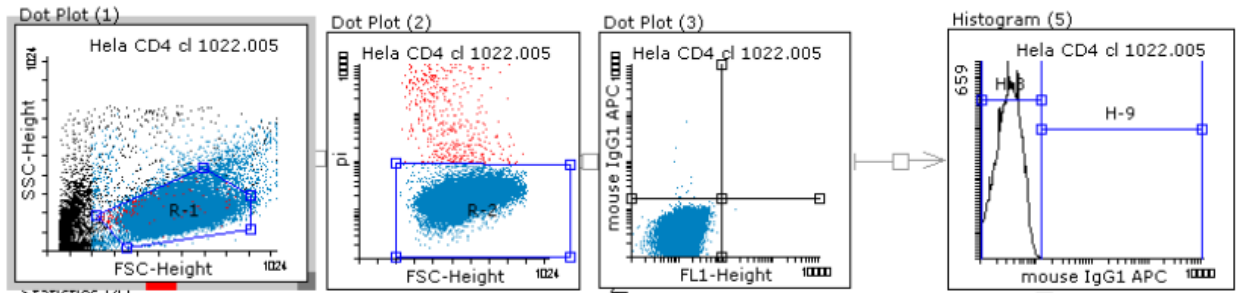


Results: 0.03% of cells showed positive fluorescence for CD4.

Figure 8: HeLa CD4+ Cells (Clone 1022) were not selected with antibiotics and CD4 was APC stained

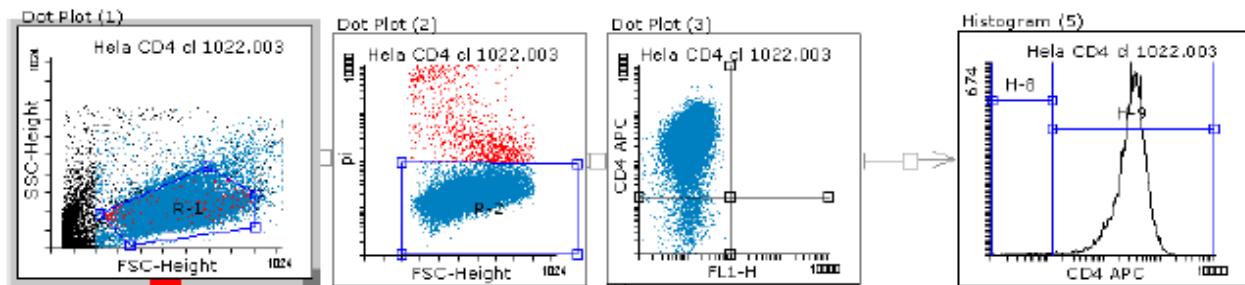


Results: 98.02% of cells showed positive fluorescence for CD4.

Figure 9: HeLa CD4+ Cells (Clone 1022) were selected with antibiotics and isotype was APC stained (Isotype Control)

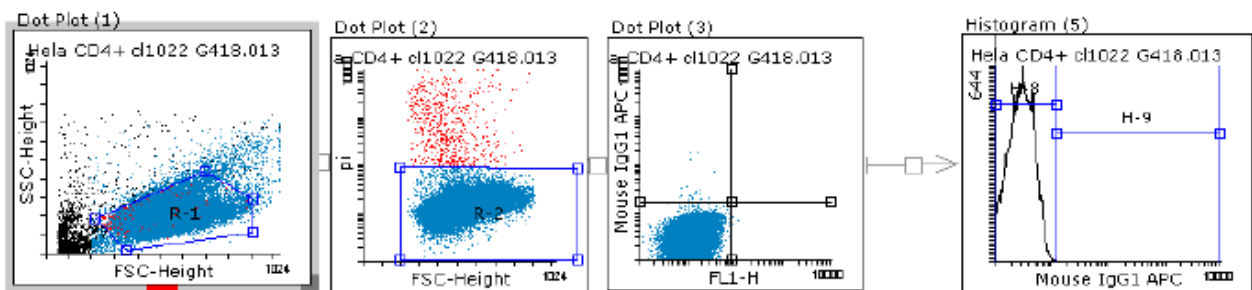

Results: 0.06% of cells showed positive fluorescence for CD4.

Figure 10: HeLa CD4+ Cells (Clone 1022) were selected with antibiotics and CD4 was APC stained

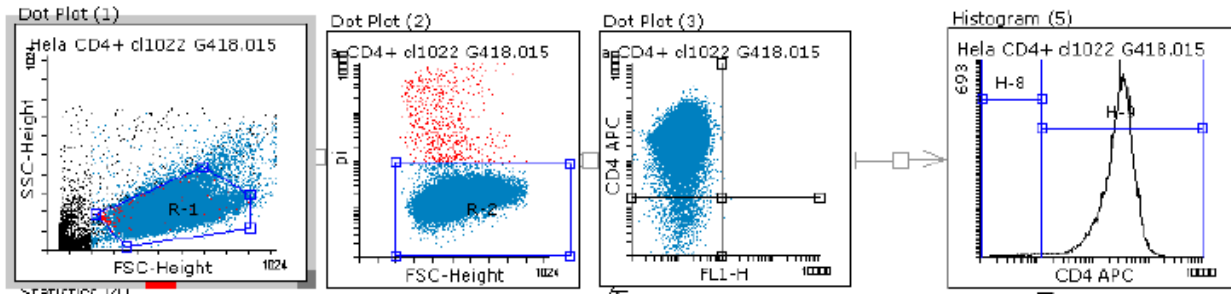

Results: 97.85% of cells showed positive fluorescence for CD4.