

Certificate of Analysis for NR-52162

Acinetobacter baumannii, Strain MRSN 3658

Catalog No. NR-52162

This reagent is the tangible property of the U.S. Government.

Product Description:

Acinetobacter baumannii (A. baumannii), strain MRSN 3658 was isolated in 2010 from a human specimen in Europe as part of a global surveillance program. A. baumannii, strain MRSN 3658 was deposited as sensitive to amikacin, ampicillin/sulbactam, cefepime, ceftazidime, ceftriaxone, ciprofloxacin, colistin, gentamicin, imipenem, levofloxacin, meropenem, tetracycline, trimethoprim/sulfamethoxazole and tobramycin. NR-52162 lot 70041722 was produced by inoculation of the deposited material into Tryptic Soy broth and grown for 1 day at 37°C in an aerobic atmosphere. Broth inoculum was added to Tryptic Soy agar kolles, which were grown for 1 day at 37°C in an aerobic atmosphere to produce this lot. Quality control testing was completed under propagation conditions unless otherwise noted.

Lot: 70041722 Manufacturing Date: 05FEB2021

BEI Resources is committed to ensuring digital accessibility for people with disabilities. This Certificate of Analysis contains complex tables and may not be fully accessible. Please let us know if you encounter accessibility barriers and a fully accessible document will be provided: E-mail: Contact@BEIResources.org. We try to respond to feedback within 24 hours.

TEST	SPECIFICATIONS	RESULTS
Phenotypic Analysis		
Cellular morphology	Gram-negative rods	Gram-negative rods
Colony morphology	Report results	Circular, convex, entire, smooth and cream (Figure 1)
Growth at 44°C ± 2°C ¹ 1 day in an aerobic atmosphere on Tryptic Soy agar	Growth	Growth
Motility	Report results	Motile
Remel™ Motility Test Medium w/TTC Indicator for 1 day at 37°C in an aerobic atmosphere		
VITEK® MS (MALDI-TOF)	A. baumannii	A. baumannii (99.9%)
Antibiotic Susceptibility Profile ^{2,3}		
Amikacin	Sensitive	Sensitive (4 to 6 µg/mL)
Ampicillin/sulbactam	Sensitive	Sensitive (1 µg/mL)
Cefepime	Sensitive	Sensitive (1.5 µg/mL)
Ceftriaxone	Sensitive	Intermediate (12 to 16 µg/mL) ⁴
Ceftazidime	Sensitive	Sensitive (4 µg/mL)
Ciprofloxacin	Sensitive	Sensitive (≤ 0.25 µg/mL)
Colistin	Sensitive	Sensitive (≤ 1 µg/mL)
Gentamicin	Sensitive	Sensitive (≤ 1 µg/mL)
Imipenem	Sensitive	Sensitive (≤ 1 µg/mL)
Levofloxacin	Sensitive	Sensitive (≤ 1 μg/mL)
Meropenem	Sensitive	Sensitive (0.38 µg/mL)
Trimethoprim/sulfamethoxazole	Sensitive	Sensitive (≤ 2 µg/mL)
Tobramycin	Sensitive	Sensitive (≤ 1 µg/mL)
Tetracycline	Sensitive	Sensitive (1.5 µg/mL)
Genotypic Analysis Sequencing of 16S ribosomal RNA gene (1470 base pairs)	≥ 99% sequence identity to A. baumannii, strain MRSN 3658 (GenBank: VHEV01000059.1)	100% sequence identity to A. baumannii, strain MRSN 3658 (GenBank: VHEV01000059.1)
Purity 7 days at 37°C in an aerobic atmosphere with and without 5% CO ₂ on Tryptic Soy agar	Growth consistent with expected colony morphology	Growth consistent with expected colony morphology

BEI Resources www.beiresources.org E-mail: contact@beiresources.org
Tel: 800-359-7370

Fax: 703-365-2898




Certificate of Analysis for NR-52162

TEST	SPECIFICATIONS	RESULTS
Viability	Growth	Growth

¹Growth at 44°C differentiates A. baumannii from A. calcoaceticus and A. pittii, which do not grow at 44°C.

Figure 1: Colony Morphology

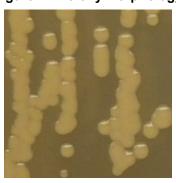

/Heather Couch/

Heather Couch 22 NOV 2021

Program Manager or designee, ATCC Federal Solutions

ATCC®, on behalf of BEI Resources, hereby represents and warrants that the material provided under this certificate has been subjected to the tests and procedures specified and that the results described, along with any other data provided in this certificate, are true and accurate to the best of ATCC®'s knowledge.

ATCC® is a trademark of the American Type Culture Collection.

You are authorized to use this product for research use only. It is not intended for human use.

BEI Resources www.beiresources.org E-mail: contact@beiresources.org
Tel: 800-359-7370

Fax: 703-365-2898

²Minimum Inhibitory Concentration (MIC); MIC Interpretation Guideline: CLSI M100-S28 (2018)

³Antibiotic susceptibility was tested using a combination of VITEK2 GN82, Sensititre GNX2F AST and E-test strips.

⁴The susceptibility result for this antibiotic is within one doubling dilution of specification, which is considered an equivalent result.