

Certificate of Analysis for NR-52188

Acinetobacter baumannii, Strain MRSN 11695

Catalog No. NR-52188

This reagent is the tangible property of the U.S. Government.

Product Description:

Acinetobacter baumannii (A. baumannii), strain MRSN 11695 was isolated in 2007 from a human urine sample in the United States as part of a global surveillance program. A. baumannii, strain MRSN 11695 was deposited as sensitive to colistin, intermediately resistant to tetracycline and resistant to amikacin, ceftazidime, ciprofloxacin, ceftriaxone, cefepime, gentamicin, imipenem, tobramycin, trimethoprim/sulfamethoxazole, levofloxacin, meropenem and ampicillin/sulbactam. NR-52188 was produced by inoculation of BEI Resources seed lot 70039039 into Tryptic Soy broth and grown for 1 day at 37°C in an aerobic atmosphere. Broth inoculum was added to Tryptic Soy agar kolles, which were grown for 1 day at 37°C in an aerobic atmosphere to produce this lot. Quality control testing was completed under propagation conditions unless otherwise noted.

Lot: 70053517 Manufacturing Date: 16JUN2022

BEI Resources is committed to ensuring digital accessibility for people with disabilities. This Certificate of Analysis contains complex tables and may not be fully accessible. Please let us know if you encounter accessibility barriers and a fully accessible document will be provided: E-mail: Contact@BEIResources.org. We try to respond to feedback within 24 hours.

TEST	SPECIFICATIONS	RESULTS
Phenotypic Analysis		
Cellular morphology	Gram-negative rods	Gram-negative rods
Colony morphology	Report results	Circular, convex, entire, smooth, mucoid and cream (Figure 1)
Growth at 44°C + 2°C ¹	Growth	Growth
1 day in an aerobic atmosphere on Tryptic	Glowali	Glowal
Soy agar		
Motility	Non-motile	Non-motile
BD BBL™ Motility Test Medium with TTC		
Indicator for 1 day at 37°C in an aerobic		
atmosphere		
VITEK® MS (MALDI-TOF)	A. baumannii	A. baumannii (99.9%)
Antibiotic Susceptibility Profile ^{2,3}		
Amikacin	Resistant	Resistant (≥ 256 µg/mL)
Ampicillin/sulbactam	Sensitive	Sensitive (8 μg/mL) ⁴
Cefepime	Resistant	Resistant (≥ 256 µg/mL)
Ceftazidime	Resistant	Resistant (≥ 64 µg/mL)
Ceftriaxone	Resistant	Resistant (≥ 64 µg/mL)
Ciprofloxacin	Resistant	Resistant (≥ 4 µg/mL)
Gentamicin	Resistant	Resistant (≥ 16 µg/mL)
Imipenem	Resistant	Resistant (≥ 32 µg/mL)
Levofloxacin	Resistant	Resistant (≥ 8 µg/mL)
Meropenem	Resistant	Resistant (≥ 16 µg/mL)
Tetracycline	Resistant	Resistant (≥ 256 µg/mL) ⁵
Tobramycin	Resistant	Resistant (≥ 16 µg/mL)
Trimethoprim/sulfamethoxazole	Resistant	Resistant (160 µg/mL)
Genotypic Analysis		
Sequencing of 16S ribosomal RNA gene	≥ 99% sequence identity to	99.9% sequence identity to
(~ 1450 base pairs)	A. baumannii, strain MRSN 11695	A. baumannii, strain MRSN 11695
	(GenBank: VHHI01000080.1)	(GenBank: VHHI01000080.1)
Purity	Growth consistent with expected	Growth consistent with expected colony
7 days at 37°C in an aerobic atmosphere with	colony morphology	morphology
and without 5% CO ₂ on Tryptic Soy agar		

BEI Resources www.beiresources.org E-mail: contact@beiresources.org

Tel: 800-359-7370 Fax: 703-365-2898




Certificate of Analysis for NR-52188

TEST	SPECIFICATIONS	RESULTS
Viability	Growth	Growth

¹Growth at 44°C differentiates A. baumannii from A. calcoaceticus and A. pittii, which do not grow at 44°C.

Figure 1: Colony Morphology

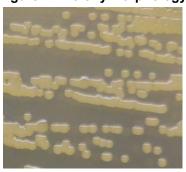

/Sonia Bjorum Brower/ Sonia Bjorum Brower

20 SEP 2023

Technical Manager or designee, ATCC Federal Solutions

ATCC®, on behalf of BEI Resources, hereby represents and warrants that the material provided under this certificate has been subjected to the tests and procedures specified and that the results described, along with any other data provided in this certificate, are true and accurate to the best of ATCC®'s knowledge.

ATCC® is a trademark of the American Type Culture Collection.

You are authorized to use this product for research use only. It is not intended for human use.

BEI Resources www.beiresources.org E-mail: contact@beiresources.org
Tel: 800-359-7370

Fax: 703-365-2898

²Minimum Inhibitory Concentration (MIC); MIC Interpretation Guideline: CLSI M100-S28 (2018)

³Antibiotic susceptibility was tested using a combination of VITEK®2 GN81 and E-test strips.

⁴A. baumannii, strain MRSN 11695 was deposited as resistant to ampicillin/sulbactam, but showed MICs of 4 to 6 μg/mL (interpreted as sensitive) for lot 70039038 during QC testing.

⁵A. baumannii, strain MRSN 11695 was deposited as intermediately resistant to tetracycline, but showed a MIC of ≥ 256 μg/mL (interpreted as resistant) for lot 70039038 during QC testing.