

Staphylococcus aureus, Strain HT 20020320

Catalog No. NR-46038

Product Description: *Staphylococcus aureus* (*S. aureus*), strain HT 20020320 was isolated in 2002 from a wound of an 8-year-old female with impetigo in France. *S. aureus*, strain HT 20020320 is a clinically associated methicillin-resistant *S. aureus* (MRSA) strain.


Lot¹: 63397520

Manufacturing Date: 27MAR2015

TEST	SPECIFICATIONS	RESULTS
Phenotypic Analysis Cellular morphology Colony morphology ² Motility (wet mount) Hemolysis ² Biochemical Analysis Catalase Coagulase ³ VITEK [®] 2 Compact (GP card)	Gram-positive cocci Report results Report results Report results Positive Report results Consistent with <i>S. aureus</i>	Gram-positive cocci Circular, low convex, entire, smooth and cream (Figure 1) Non-motile β-hemolytic Positive Positive Consistent with <i>S. aureus</i>
Antibiotic Susceptibility Profile VITEK [®] (AST-GP71 card) ⁴ Beta-lactamase ⁵ Cefoxitin screen Benzylpenicillin Oxacillin Gentamicin Ciprofloxacin Levofloxacin Moxifloxacin Clindamycin (inducible resistance) Erythromycin Clindamycin Quinupristin/dalfopristin Linezolid Daptomycin Vancomycin Minocycline Tetracycline Tigecycline Nitrofurantoin Rifampicin Trimethoprim/sulfamethoxazole Etest [®] antibiotic test strips ⁶ Chloramphenicol ⁷ Teicoplanin ⁷	Report results Report results Report results Resistant Sensitive Resistant Report results Report results Report results Sensitive Sensitive Sensitive Sensitive Report results Sensitive Report results Report results Report results Report results Sensitive Report results Sensitive	Positive Positive Resistant (≥ 0.5 µg/mL) Resistant (≥ 4 µg/mL) Sensitive (≤ 0.5 µg/mL) Resistant (≥ 8 µg/mL) Resistant (= 4 µg/mL) Intermediate (= 1 µg/mL) Negative Sensitive (≤ 0.25 µg/mL) Sensitive (≤ 0.25 µg/mL) Sensitive (≤ 0.25 µg/mL) Sensitive (= 2 µg/mL) Sensitive (= 0.5 µg/mL) Sensitive (= 1 µg/mL) Sensitive (≤ 0.5 µg/mL) Sensitive (≤ 1 µg/mL) Sensitive (≤ 0.12 µg/mL) Sensitive (≤ 16 µg/mL) Sensitive (≤ 0.5 µg/mL) Sensitive (≤ 10 µg/mL) Sensitive (= 3-4 µg/mL) Sensitive (= 1.5-2 µg/mL)
Genotypic Analysis Sequencing of 16S ribosomal RNA gene (~ 1500 base pairs)	Consistent with <i>S. aureus</i>	Consistent with <i>S. aureus</i>
Purity (post-freeze)⁸	Growth consistent with <i>S. aureus</i>	Growth consistent with <i>S. aureus</i>
Viability (post-freeze)²	Growth	Growth

- ¹*S. aureus*, strain HT 20020320 was deposited to BEI Resources as part of the NARSA collection. NR-46038 was produced by inoculation of the deposited material into Tryptic Soy broth and grown 24 hours at 37°C in an aerobic atmosphere. Broth inoculum was added to Tryptic Soy agar with 5% defibrinated sheep blood kolles which were grown 24 hours at 37°C in an aerobic atmosphere to produce this lot.
- ²19 hours at 37°C in an aerobic atmosphere on Tryptic Soy agar with 5% defibrinated sheep blood
- ³4 hours at 37°C in rabbit serum with 0.15% EDTA (Coagulase Plasma BBL™ 240827)
- ⁴Minimum Inhibitory Concentration (MIC); MIC Interpretation Guideline: CLSI M100-S22 (2012)
- ⁵The production of beta-lactamase was detected using a Cefinase™ Paper Disc (BBL™ 231650).
- ⁶24 hours at 37°C in an aerobic atmosphere on Mueller Hinton agar
- ⁷For both chloramphenicol (bioMérieux Etest® 412308) and teicoplanin (bioMérieux Etest® 412459), a MIC ≤ 8 µg/mL is sensitive, a MIC = 16 µg/mL is intermediate and a MIC ≥ 32 µg/mL is resistant.
- ⁸Purity of this lot was assessed for 7 days on Tryptic Soy agar with 5% defibrinated sheep blood at 37°C in an aerobic atmosphere.

Figure 1

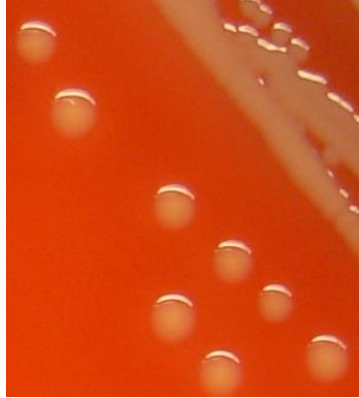

Date: 27 MAY 2015

Signature: 

BEI Resources Authentication

ATCC®, on behalf of BEI Resources, hereby represents and warrants that the material provided under this certificate has been subjected to the tests and procedures specified and that the results described, along with any other data provided in this certificate, are true and accurate to the best of ATCC®'s knowledge.

ATCC® is a trademark of the American Type Culture Collection.

You are authorized to use this product for research use only. It is not intended for human use.

