

Certificate of Analysis for NR-55506

Klebsiella pneumoniae, Strain MRSN 4759

Catalog No. NR-55506

This reagent is the tangible property of the U.S. Government.

Product Description:

Klebsiella pneumoniae (K. pneumoniae), strain MRSN 4759 was isolated in 2011 from a human urine sample in North America as part of a global surveillance program. NR-55506 was deposited as a multidrug-resistant strain, sensitive to amikacin, cefepime, ceftazidime/avibactam, ceftolozane/tazobactam, ciprofloxacin, ertapenem, gentamicin, imipenem, levofloxacin, meropenem, piperacillin/tazobactam, tetracycline, tigecycline, tobramycin and trimethoprim/sulfamethoxazole and resistant to ampicillin/sulbactam, aztreonam, ceftazidime and ceftriaxone. NR-55506 was produced by inoculation of the deposited material into Tryptic Soy broth and grown for 1 day at 37°C in an aerobic atmosphere. Broth inoculum was added to Tryptic Soy agar kolles, which were grown for 1 day at 37°C in an aerobic atmosphere to produce this lot. Quality control testing was completed under propagation conditions unless otherwise noted.

Lot: 70049672 Manufacturing Date: 20JAN2022

BEI Resources is committed to ensuring digital accessibility for people with disabilities. This Certificate of Analysis contains complex tables and may not be fully accessible. Please let us know if you encounter accessibility barriers and a fully accessible document will be provided: E-mail: contact@BEIResources.org. We try to respond to feedback within 24 hours.

TEST	SPECIFICATIONS	RESULTS
Phenotypic Analysis		
Cellular morphology	Gram-negative rods	Gram-negative rods
Colony morphology	Report results	Circular, convex, entire, smooth, mucoid and cream (Figure 1)
Motility (wet mount)	Report results	Non-motile
VITEK® 2 (GN card)	K. pneumoniae (≥ 89%)	K. pneumoniae (99%)
VITEK 2 (GN card) VITEK® MS (MALDI-TOF)	K. pneumoniae	K. pneumoniae (99.9%)
Antibiotic Susceptibility Profile ^{1,2}	N. prieumomae	N. prieumoniae (99.970)
Amikacin	Sensitive	Sensitive (≤ 2 µg/mL)
Ampicillin/sulbactam	Resistant	Resistant (≥ 32 µg/mL)
Aztreonam	Resistant	Sensitive (2 µg/mL) ³
Cefepime	Sensitive	Sensitive (2 µg/mL)
Ceftazidime	Resistant	Sensitive (2 µg/mL) ⁴
Ceftazidime/avibactam	Sensitive	Sensitive (1 µg/mL)
Ceftolozane/tazobactam	Sensitive	Sensitive (0.19 µg/mL)
Ceftriaxone	Resistant	Resistant (≥ 64 µg/mL)
	Sensitive	Sensitive (1 µg/mL)
Ciprofloxacin	Sensitive	Sensitive (1 µg/mL) Sensitive (≤ 0.5 µg/mL)
Ertapenem Gentamicin	Sensitive	Sensitive (≤ 0.5 µg/mL) Sensitive (≤ 1 µg/mL)
	Sensitive	` ,
lmipenem Levofloxacin	Sensitive	Sensitive (0.25 µg/mL)
		Sensitive (1 µg/mL)
Meropenem	Sensitive	Sensitive (≤ 0.25 μg/mL)
Piperacillin/tazobactam	Sensitive	Sensitive (≤ 4 μg/mL)
Tetracycline	Sensitive	Sensitive (≤ 1 µg/mL)
Tigecycline	Sensitive	Sensitive (≤ 0.5 µg/mL) ⁵
Tobramycin	Sensitive	Sensitive (≤ 1 μg/mL)
Trimethoprim/sulfamethoxazole	Sensitive	Sensitive (≤ 20 μg/mL)
Genotypic Analysis	> 000/ apquence identity to	00.3% acquence identity to
Sequencing of 16S ribosomal RNA gene (~ 1470 base pairs)	≥ 99% sequence identity to K. pneumoniae, strain MRSN 4759 (GenBank: JAGYFA010000074.1)	99.3% sequence identity to K. pneumoniae, strain MRSN 4759 (GenBank: JAGYFA010000074.1) ⁶

BEI Resources www.beiresources.org E-mail: contact@beiresources.org

Tel: 800-359-7370 Fax: 703-365-2898




Certificate of Analysis for NR-55506

TEST	SPECIFICATIONS	RESULTS
Purity 7 days at 37°C in an aerobic atmosphere with 5% CO ₂ on Tryptic Soy agar with 5% defibrinated sheep blood	colony morphology	Growth consistent with expected colony morphology
Viability	Growth	Growth

¹Minimum Inhibitory Concentration (MIC); MIC Interpretation Guideline: CLSI M100-S28 (2018)

Figure 1: Colony Morphology

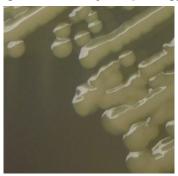

/Heather Couch/

Heather Couch 18 APR 2022

Program Manager or designee, ATCC Federal Solutions

ATCC®, on behalf of BEI Resources, hereby represents and warrants that the material provided under this certificate has been subjected to the tests and procedures specified and that the results described, along with any other data provided in this certificate, are true and accurate to the best of ATCC®'s knowledge.

ATCC® is a trademark of the American Type Culture Collection.

You are authorized to use this product for research use only. It is not intended for human use.

BEI Resources www.beiresources.org E-mail: contact@beiresources.org Tel: 800-359-7370

Fax: 703-365-2898

²Antibiotic susceptibility was tested using a combination of bioMérieux VITEK®2 GN74 and ETÉST®.

³K. pneumoniae, strain MRSN 4759 was deposited as resistant to aztreonam, but showed a MIC of 2 µg per mL (interpreted as sensitive) for this antibiotic during QC testing. Testing was performed in duplicate.

4K. pneumoniae, strain MRSN 4759 was deposited as resistant to ceftazidime, but showed a MIC of 1 µg per mL (interpreted as sensitive) for this

antibiotic during QC testing. Testing was performed in duplicate.

MIC Interpretation Guideline: EUCAST Version 8.0 (2018)

⁶Also consistent with other Klebsiella species