

***Mycobacterium tuberculosis*, Strain KT-0037**

**Catalog No. NR-43822**

**Product Description:**

*Mycobacterium tuberculosis* (*M. tuberculosis*), strain KT-0037 was isolated from a human in South Korea. Strain KT-0037 was deposited as a multidrug-resistant (MDR) Beijing genotype strain with resistance to resistance to capreomycin, isoniazid, kanamycin, moxifloxacin, ofloxacin, pyrazinamide, rifampin, pyrazinamide and streptomycin. NR-43822 was produced by inoculation of the deposited material into Middlebrook 7H9 broth with ADC enrichment. Broth inoculum was added to Middlebrook 7H10 agar with OADC enrichment kolles, which were grown for 21 days at 37°C in an aerobic atmosphere with 5% CO<sub>2</sub> to produce this lot.

**Lot: 70023004**

**Manufacturing Date: 03APR2019**

BEI Resources is committed to ensuring digital accessibility for people with disabilities. This Certificate of Analysis contains complex tables and may not be fully accessible. Please let us know if you encounter accessibility barriers and a fully accessible document will be provided: E-mail: [Contact@BEIResources.org](mailto:Contact@BEIResources.org). We try to respond to feedback within 24 hours.

TEST	SPECIFICATIONS	RESULTS
<b>Phenotypic Analysis<sup>1</sup></b> Cellular morphology 21 days at 37°C in an aerobic atmosphere with 5% CO <sub>2</sub> on Middlebrook 7H10 agar with OADC enrichment Colony morphology Growth rate Growth at 26°C Growth at 37°C Acid-fast stain Pigmentation in the dark (Scotochromogen) Photoinduction for 1 hour (Photochromogen) Nonchromogen (no pigment) Biochemical tests Niacin production <sup>2</sup> Nitrate reduction	Gram-positive rods Report results ≥ 7 days Negative Positive Positive (red colonies) Negative (no pigment) Negative (no pigment) Positive (no pigment) Positive Positive	Gram-positive rods Irregular, slight peaked, undulate, rough and cream (Figure 1) 21 days Negative Positive Positive (red colonies) Negative (no pigment) Negative (no pigment) Positive (no pigment) Positive Positive
<b>Antibiotic Susceptibility Profile</b> Sensititre™ System <sup>3,4</sup> Amikacin Cycloserine Ethambutol Ethionamide Isoniazid Kanamycin Moxifloxacin Ofloxacin Para-aminosalicylic acid Rifabutin Rifampin Streptomycin	Report results Report results Report results Report results Report results Report results Report results Report results Report results Report results Report results	16 µg/mL 32 µg/mL <sup>5,6</sup> 8 µg/mL <sup>7</sup> 20 µg/mL <sup>6,7,8</sup> > 4 µg/mL ≥ 40 µg/mL 4 µg/mL <sup>6,9</sup> 8 µg/mL > 64 µg/mL <sup>7</sup> 8 µg/mL <sup>6,7,10</sup> > 16 µg/mL 2 µg/mL <sup>6,7,11</sup>
<b>Genotypic Analysis</b> Sequencing of Heat Shock Protein 65 gene (~ 1620 base pairs)	≥ 99% sequence identity to <i>M. tuberculosis</i> , strain KT-0037 (GenBank: JLN101000007.1)	100% sequence identity to <i>M. tuberculosis</i> , strain KT-0037 (GenBank: JLN101000007.1 <sup>12</sup> )

TEST	SPECIFICATIONS	RESULTS
<b>Purity (post-freeze)</b> Middlebrook 7H10 agar with OADC enrichment 37 days at 37°C in an aerobic atmosphere with 5% CO <sub>2</sub> Tryptic Soy agar 21 days at 37°C in an aerobic atmosphere with 5% CO <sub>2</sub>	Growth consistent with expected colony morphology  Report results	Growth consistent with expected colony morphology  Growth consistent with expected colony morphology
<b>Viability (post-freeze)</b> 21 days at 37°C in an aerobic atmosphere with 5% CO <sub>2</sub> on Middlebrook 7H10 agar with OADC enrichment	Growth	Growth

<sup>1</sup>Information on *Mycobacterium* testing is available from Ribón, W. "Biochemical Isolation and Identification of Mycobacteria." Biochemical Testing. (2012) Jose C. Jimenez-Lopez (Ed.), InTech ([Biochemical Isolation and Identification of Mycobacteria](#)) and Lévy-Frébault, V. V. and F. Portaels. "Proposed Minimal Standards for the Genus *Mycobacterium* and for Description of New Slowly Growing *Mycobacterium* Species." Int. J. Syst. Bacteriol. 42 (1992): 315-323. PubMed: 1581193.

<sup>2</sup>All mycobacteria produce niacin but only *M. tuberculosis* accumulates it, resulting in a positive test for *M. tuberculosis*.

<sup>3</sup>Sensititre™ System *Mycobacterium tuberculosis* MIC Plate, Thermo Scientific™, catalog number MYCOTB

<sup>4</sup>Minimum Inhibitory Concentration (MIC); No Clinical & Laboratory Standards Institute (CLSI) interpretations of the Sensititre™ System data for *M. tuberculosis* are currently available.

<sup>5</sup>Two MICs were observed for cycloserine (16 µg/mL and 32 µg/mL) under identical test conditions. The highest MIC is being reported as the test result.

<sup>6</sup>Variability in the MIC result by the Sensititre™ method has been demonstrated [Lee, J., et al. "Sensititre MYCOTB MIC Plate for Testing *Mycobacterium tuberculosis* Susceptibility to First- and Second-Line Drugs." Antimicrob. Agents Chemother. 58 (2014): 11-18. PubMed: 24100497.], with the results for a single antibiotic typically within one doubling dilution.

<sup>7</sup>For ethambutol, ethionamide, para-aminosalicylic acid, rifabutin and streptomycin, the endpoint for these drugs is determined by the well with approximately 80% inhibition of growth compared to the positive control well with no drug.

<sup>8</sup>Two MICs were observed for ethionamide (10 µg/mL and 20 µg/mL) under identical test conditions. The highest MIC is being reported as the test result.


<sup>9</sup>Two MICs were observed for moxifloxacin (2 µg/mL and 4 µg/mL) under identical test conditions. The highest MIC is being reported as the test result.

<sup>10</sup>Two MICs were observed for rifabutin (4 µg/mL and 8 µg/mL) under identical test conditions. The highest MIC is being reported as the test result.

<sup>11</sup>Two MICs were observed for streptomycin (1 µg/mL and 2 µg/mL) under identical test conditions. The highest MIC is being reported as the test result.

<sup>12</sup>Also consistent with *M. africanum*, *M. bovis*, *M. canettii*, *M. caprae* and *M. microti*

Figure 1: Colony Morphology

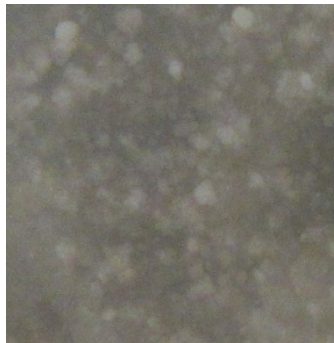

/Heather Couch/  
 Heather Couch

30 MAR 2022

Program Manager or designee, ATCC Federal Solutions

ATCC®, on behalf of BEI Resources, hereby represents and warrants that the material provided under this certificate has been subjected to the tests and procedures specified and that the results described, along with any other data provided in this certificate, are true and accurate to the best of ATCC®'s knowledge.

ATCC® is a trademark of the American Type Culture Collection.  
 You are authorized to use this product for research use only. It is not intended for human use.

